# BRAUN

## Oral-B

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## English

Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

### IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the toothbrush, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

### WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

### Description

- A Brush head (depending on model/ interchangeable)
- B On/off switch
- C Handle
- D Charging unit

### Specifications

For voltage specifications please refer to the bottom of the charging unit.

### Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the handle (C) on there (1). For voltage specifications refer to the bottom of the charging unit.
- A full charge can take up to 16 hours and enables up to 5 days of regular brushing (twice a day, 2 minutes).

- For everyday use, the handle can be stored on the plugged-in charging unit. There is no risk of overcharging.

### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

### Using the toothbrush

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (2).

When using one of the rotating brush heads guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface (3).

When using the TriZone brush head place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush.

Start brushing the outsides, then the insides, finally the chewing surfaces. Brush at least 2 minutes to thoroughly remove plaque. Press the on/off switch to turn off.

### Individual brush heads

Oral-B offers you a variety of brush heads which are available in stores that sell Oral-B products. Below brush heads are interchangeable and compatible with your Oral-B handle.

#### The Oral-B Precision Clean brush head

cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

#### The Oral-B FlossAction brush head

effectively cleans your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth.

#### The Oral-B ProWhite brush head

features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening.

Please note that children under 12 years old should not use the Oral-B ProWhite brush head.

Our dynamic TriZone brush head can be used with all common brushing techniques.

#### The Oral-B TriZone brush head

features triple zone cleaning action for outstanding plaque removal, even in between teeth.

We do not recommend using the Oral-B FlossAction or Oral-B ProWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

All brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need (3.a). With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months (3.b), signaling it is time to change the brush head. If the bristles splay before the colour recedes (3.c), you are exerting too much pressure on teeth and gums.

During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.

### Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on. Switch off and remove the brush head. Rinse both parts separately under running water and wipe them dry. From time to time, also clean the charging unit with a damp cloth (4).

Subject to change without notice.

### Environmental notice


This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at

the end of its useful life. Disposal can take place at a Braun Oral-B Service Centre or at appropriate collection points provided in your country.

However, if you prefer to dispose of the battery yourself, open the handle as shown below:

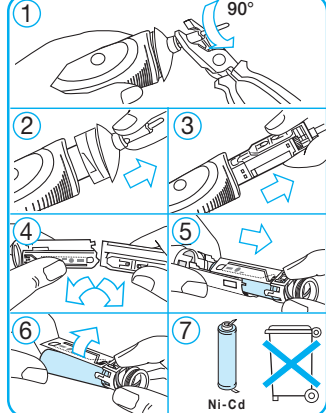**Caution! Opening the handle will destroy the appliance and invalidate the guarantee.**

### Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase (Guarantee period may differ in different country). Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Oral-B Service Centre.

**Charger type 3757**  
**Handle type 3709**

96186133/IV-13  
Japan/GB