BRAUN

Oral-B

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Braun Oral-B Denta-Pride has been carefully designed to offer you and your family a new and unique brushing experience that is both safe and effective. However, when using electrical products, basic safety precautions should always be followed. Please read carefully before use.

Important

- **Braun Oral-B Denta-Pride** is provided with a SmartPlug with integrated Safety Extra Low Voltage power supply. Do not exchange or manipulate any part of it. Otherwise there is the risk of electric shock. Periodically check the SmartPlug for damage. A damaged or non-functioning appliance shall no longer be used. Take it to a Braun Oral-B Repair Centre.
- This product is not intended for use by children under age three.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.

Description of Braun Oral-B Denta-Pride

- a Braun Oral-B FlossAction brush head
- b On/off button
- c Brushing mode button
- d Handle
- e SmartDisplay
- f SmartPlug charger (includes charging ring and SmartPlug)
- g Base plate
- h Base station
- i Protective cover
- j Braun Oral-B Professional White brush head (deluxe model only)
- k Tongue Freshener


A convenient travel case (deluxe model only) holds the **Braun Oral-B Denta-Pride** handle and two Braun Oral-B brush heads (except for Braun Oral-B CrossAction Power and Braun Oral-B Sonic Complete brush heads).

Connecting and charging

Braun Oral-B Denta-Pride has a waterproof handle, is electrically safe and is designed for use in the bathroom.

To assemble the base station (h) as a charging unit, detach the base plate (g) first (see diagram ①, page 3). Place the charging ring into the base station and attach the base plate as shown (② and ③, page 3). Then insert the protective cover (i) on the back of your base station (④, page 3).


- Plug the SmartPlug into an electrical outlet and place the toothbrush handle (d) on the charging unit. Alternatively, **Braun Oral-B Denta-Pride** can be charged by placing the handle directly into the charging ring while the SmartPlug is plugged in.

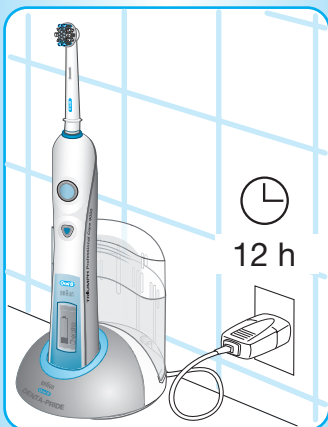- A battery symbol should appear in the SmartDisplay. The battery symbol will fill up while the handle is being charged and «Charging» will be displayed. When the handle is fully charged, the SmartDisplay shows «Full Charge». A full charge takes at least 12 hours and enables approximately two weeks of regular brushing.

Please note: If the battery symbol is not displayed immediately, continue charging. It will appear in 10 – 15 minutes.

- When the rechargeable battery is low, the bars of the battery symbol will disappear and «Low» will be displayed. Your toothbrush will slow down and then stop. It is time to recharge **Braun Oral-B Denta-Pride**. When the battery is fully discharged, «Empty» is displayed. It will then need a 15-minute charge before you can use it for a 2-minute brushing (twice a day, 2 minutes).
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months. **Braun Oral-B Denta-Pride** is convenient for travelling because you only need the compact SmartPlug charger and an appropriate travel adapter for recharging in foreign countries. For electric specifications, see printing on the SmartPlug.



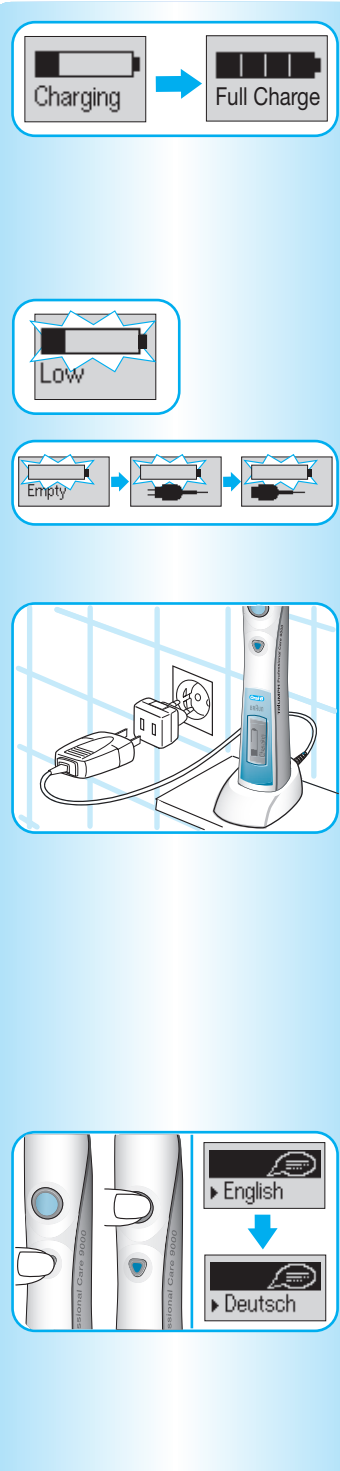
Personalizing Braun Oral-B Denta-Pride

Before getting started, you can program your toothbrush so that it meets your specific needs. Braun Oral-B Denta-Pride has innovative Smart Technology that lets you choose the language and the timer that works best for you.

Select your language

Braun Oral-B Denta-Pride comes with the English language activated. If you wish to switch to another language, proceed as follows: Press and hold the brushing mode button (c) until «English» is displayed. By briefly pressing the on/off button (b) you can scroll through 13 different languages. To confirm your desired language, briefly press the brushing mode button again. This will also lead you directly to the timer choice menu.

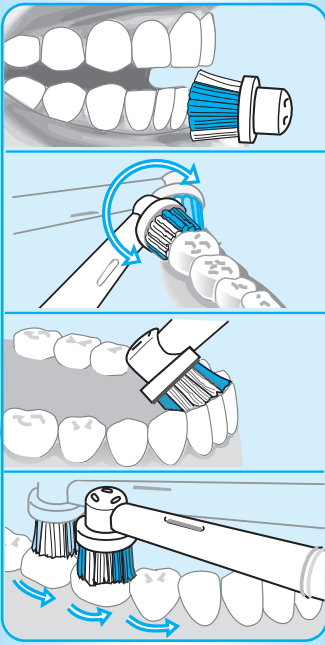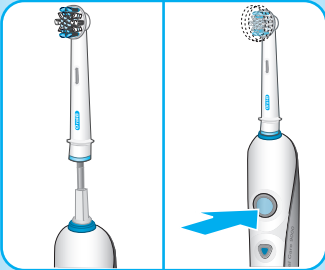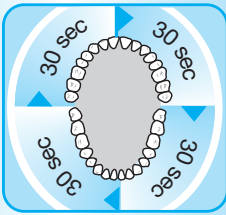Select your timer

The «2-Minute Timer» indicates by a long stuttering sound when the dentist recommended 2-minute brushing time has been reached. The «Professional Timer» helps you brush all four quadrants of your mouth equally. It signals with a short stuttering sound at 30-second intervals that it is time to move to the next quadrant of your mouth, followed by a long stuttering sound when you reach the recommended two minute mark.

Braun Oral-B Denta-Pride comes with the «Professional Timer» activated. If you wish to switch to the «2 Minute Timer» after you have confirmed the language, simply press the on/off button to change the setting. Then press the brushing mode button to confirm your choice

Using Braun Oral-B Denta-Pride

Attach a Braun Oral-B brush head to the handle. You can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on.

Brushing technique

- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

To help you get optimum brushing results, **Braun Oral-B Denta-Pride** has a pressure control feature to ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling the change, you can also hear a different sound.

Brushing modes

Braun Oral-B Denta-Pride offers different brushing modes for varying oral care needs:

- «Clean» – Exceptional mouth cleaning.
- «Soft» – Gentle, yet thorough cleaning for sensitive areas.
- «Massage» – Gentle stimulation of gums.
- «Polish» – Exceptional polishing and natural whitening for occasional or everyday use (deluxe model only).

How to switch between brushing modes:

- **Braun Oral-B Denta-Pride** automatically starts in the «Clean» mode when pressing the on/off button.
- To shift from the «Clean» to the «Soft» mode, press the brushing mode button once.
- If you wish to switch to the «Massage» mode, press the brushing mode button again and let the brush head pulse against your gum line.
- The «Polish» mode (deluxe model only) is best used with a Professional White brush head. Once the Professional White brush head is attached, press the on/off button and your toothbrush automatically starts in the «Polish» mode (please see «Brush heads»).
- To return from any mode to the «Clean» mode, press and hold the brushing mode button.
- To turn your brush off, press the on/off button.



During brushing, the SmartDisplay shows the brushing mode you have selected, followed by your actual brushing time. When you finish brushing, the total brushing time will be displayed for 30 seconds before the display switches back to the battery charge level. To directly get to the battery charge level display, you can clear the brushing time display by pushing the brushing mode button. In the first few days of using **Braun Oral-B Denta-Pride**, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Soft» mode.


Brush heads & accessories

The Braun Oral-B FlossAction and Professional White brush heads contain a «smart chip» to communicate with the electronics inside the **Braun Oral-B Denta-Pride** handle. The «smart chip» works for you, automatically switching to the appropriate brushing mode for the attached brush head. Nevertheless, you can select other brushing modes to customize your routine.

Braun Oral-B FlossAction brush head

The Braun Oral-B FlossAction brush head is recommended for everyday deep cleaning. It features Micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing plaque. The brush head's outer bristles surround each tooth, providing exceptional overall plaque removal.

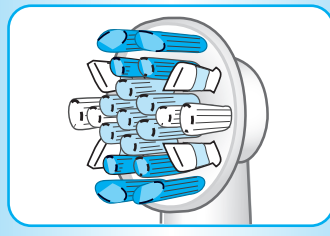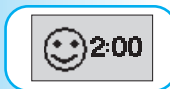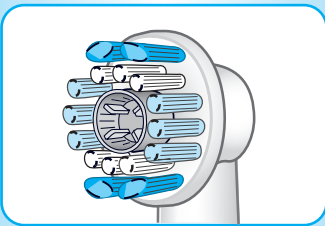Replace
Brushhead

Braun Oral-B Professional White brush head (deluxe model only)
The Braun Oral-B Professional White brush head is for occasional or everyday polishing. It features a dentist-inspired polishing cup to remove stains from coffee, tea and tobacco for effective brightening.

We do not recommend using the Braun Oral-B Professional White or Braun Oral-B FlossAction brush head with braces. The Braun Oral-B Ortho brush head is specifically designed to clean around braces and wires and can be used with the **Braun Oral-B Denta-Pride** toothbrush.

Replacement brush heads



Dentists recommend that you change your brush head every three months. As soon as your Braun Oral-B FlossAction or Professional White brush head has surpassed the recommended usage, the SmartDisplay shows «Replace Brushhead». Replace with a new brush head to get optimum brushing results. If you notice unusual wear before receiving the «Replace Brush head» message, replace your brush head immediately.

Tongue Freshener

The Tongue Freshener helps freshen breath by removing odor-causing bacteria from the tongue. It can be used everyday – before or after brushing. You can use mouthwash with the Tongue Freshener. Attach it to the handle, turn the appliance on and drag it across the tongue back to front. We recommend a low speed setting for enhanced control during use. Rinse after use. Replace the Tongue Freshener when materials become worn.


Braun Oral-B Denta-Pride replacement brush heads and accessories are available in stores that sell Braun Oral-B products.

Cleaning recommendations

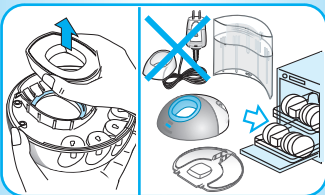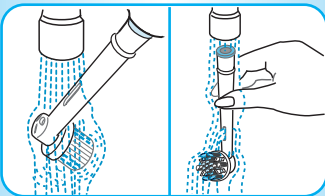- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Then wipe them dry before you put the handle on the charging unit.
- Before cleaning your charging unit, make sure you disassemble it first. The charging ring and the SmartPlug should never be placed in the dishwasher. However, all other parts (base station, base plate and protective cover) are dishwasher-safe.

- From time to time, you can also clean the charging ring and the SmartPlug with a damp cloth. The SmartPlug charger should never be placed in water.

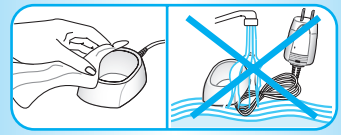Environmental notice


This product is provided with a nickel-hydride battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Oral-B Repair Centre or at approved recycling or disposing locations in accordance with local regulations.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the guarantee.

Subject to change without notice.

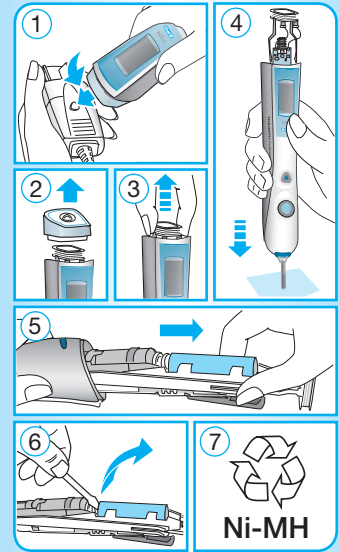Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

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