**Oral-B**

**iO**

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# Description

- 1 Brush head
- 2 Brush head identification symbol
- 3 Smart Pressure Sensor / Light Ring
- 4 Interface
- 5 *Bluetooth®*
- 6 Handle
- 7 Charger
- 8 On/off button
- 9 Mode switch button
- 10 Low charge indicator
- 11 Brush Replacement Reminder

## Features & Accessories (depending on model):

- a Smart Pressure Sensor
- b 2-minutes brushing reward illumination
- c Brush head holder with cap
- d Travel case
- e Brush head holder stand
- f Brush head «Ultimate Clean»
- g Brush head «Gentle Care»
- h Brush head «Radiant White»

## Oral-B App:

- i Download Oral-B App and connect to *Bluetooth®*  
(see back page)

**Note:** Content may vary based on model purchased.

# English

## Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

### IMPORTANT

- Periodically check the entire product/cord/accessories for damage. A damaged or non-functioning unit should no longer be used. If the product/accessory is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Never insert any object into any opening of the appliance.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Use this product only for its intended use as described in this manual. Do not

use attachments which are not recommended by the manufacturer.

- Only use the charger provided with your appliance.

### **WARNING**

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Avoid sharp bending, jamming or squeezing of the cord.
- This appliance contains batteries that are non-replaceable. The appliance has no user serviceable parts. Do not open and disassemble device except when disposing of the battery. Follow handling and disposal instructions of the battery/device described in this manual.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution. For hygienic reasons we recommend the use of individual brush handles and brush heads for each family member.
- Small parts may come off, keep out of reach of children.

To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Replace brush head every 3 months or sooner if brush head becomes worn.
- Clean brush head properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

## **Important Information**

- This toothbrush and pouch (accessory) contains NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your Oral-B iO toothbrush.
- Your toothbrush can be used with your smartphone (find details under «Connecting your toothbrush to your Smartphone»). To avoid electromagnetic interference and/or

compatibility conflicts, disable *Bluetooth*® on your toothbrush handle (6) before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals.

- Persons with pacemakers should always keep the toothbrush more than 15 centimeters (six inches) from the pacemaker while turned on. Anytime you suspect interference is taking place disable *Bluetooth*®.

## Specifications

For voltage specifications refer to the bottom of the charging unit.

Noise level: 64 dB (A)

# Using your toothbrush

## Charging and Operating

- To charge the battery place the toothbrush on the plugged-in charger (7).
- During charging the on/off button is blinking white. Once it is fully charged the indicator light turns off. A full charge typically takes 16 hours.
- If the battery is running low the low charge indicator (10) blinks in red. Once the battery is empty, the motor stops. To charge the battery please put the brush on the charger.
- The toothbrush can be stored on the plugged-in charger to maintain it at full charge; a battery overcharge is prevented by the device. However, for environmental purposes we recommend to keep charger in the power socket during charging only.
- The toothbrush switches off when exceeding the outside operating limits and/or standard battery current. To enable operation the toothbrush needs to be put on the charger.

**Caution:** This toothbrush contains a Li-Ion battery. Do not expose toothbrush to temperatures higher than 50 °C / 122 °F.

Note: Store the toothbrush at room temperature for optimal battery maintenance.

## Start Brushing

Place the brush head (1) on the handle (6). A small gap should remain between handle and brush head.

Wet brush head and apply toothpaste. To avoid splashing, put the brush head into your mouth before switching the toothbrush on by pressing the on/off button(8).

Guide the brush head from tooth to tooth. Brush the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head (g)). Availability of modes may vary depending on model purchased.

### Timer

The timer signals with a short stuttering sound in 30-seconds intervals to move to the next quadrant of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

The timer memorizes the elapsed brushing time, even when the toothbrush is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the toothbrush on the plugged-in charger.

### Brushing reward illumination

The completion of a 2-minutes brushing event will be rewarded with a colourful light ring illumination. (b)

Note: You may also customize the timer settings via Oral-B App.

### **Smart Pressure Sensor/ Light Ring (a)**

Your toothbrush is equipped with a Smart Pressure Sensor (3) that has various functionalities indicated via different colours of light (a):

White = Low pressure indicator / default setting

Green = optimal pressure indicator

Red = high pressure indicator

Brushing with the recommended pressure is indicated via a green light. If too much pressure is applied, a red light reminds you to reduce pressure. If the Smart Pressure Sensor illuminates in white (or your pre-selected colour), we recommend applying more pressure.

When high pressure is applied the brush head movement will be reduced (in «Daily Clean», «Intense» and «Sensitive» modes). Periodically check the operation of the Smart Pressure Sensor by pressing moderately on the brush head during use.

Note: You can replace the white default light ring colour with a choice of individual colours via the Oral-B App.

## Brushing modes (depending on model)

Your toothbrush offers different brushing modes:

«**Daily Clean**» – Standard mode for daily cleaning

«**Intense**» – For an extraordinary clean feeling (operated with a higher speed)

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Super Sensitive**» – Super gentle cleaning for extra sensitive areas

«**Whiten**» – Polishing for occasional or everyday use

«**Tongue Clean**» – Tongue cleaning for occasional or everyday use

«**Gum Care**» – Gentle massage of gums

Note: The sequence and availability of brushing modes may vary depending on model purchased.

### Mode setting:

For first use the toothbrush starts with the first mode illustrated on the package of your purchased model. To switch to other modes, successively press the mode button (9). You may change the modes during the brushing event. The toothbrush automatically starts in the mode selected last.

When operating the «Tongue Clean» mode, we recommend using a «Gentle Care» brush head (g). You can brush your tongue either with or without toothpaste. Brush the whole area of your tongue systematically in gentle movements. The recommended brushing time is 30 seconds; the motor stops after 30 seconds.

Note: You may also customize mode settings via Oral-B App.

## Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads (f/g/h) that only fit your Oral-B iO toothbrush handle. These unique brush heads are dentist-inspired to leave your mouth feeling clean and refreshed.

### When to exchange brush head

Some iO Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the bristle colour will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the colour recedes, you may exert too much pressure on teeth and gums. In general, we recommend replacing the brush head every 3 months.

## Brush Replacement Reminder

The toothbrush is equipped with a «Brush Replacement Reminder» indicating the need to exchange the brush head approximately after 3 months when used twice per day. A yellow\* light - accompanied by a vibration signal - will illuminate in the mode button (11).

To reset the «Brush Replacement Reminder» hold the mode button pressed for 2 seconds. A vibration and illumination signal confirms the reset. If no action is taken, an automatic reset will happen after several days.

\*may appear orange

Note: You may also customize the brush replacement settings via the Oral-B App.

## Connecting your toothbrush to your Smartphone/ Oral-B App

The Oral-B App offers various features:

- **Personalize your toothbrush settings** (e.g. choose your preferred light ring colour, adjust stutter timer, mode sequence or brush replacement reminder settings)
- **Track your brushing behavior** (e.g. brushing time and applied pressure)

The Oral-B App is available for smartphones running iOS or Android™ operating systems. It can be downloaded for free from App Store<sup>(SM)</sup>/ Google Play™ or you may scan the QR code on the back page of this manual (i).

- Launch the Oral-B App. It will guide you through the entire *Bluetooth*® pairing procedure. You may also find it under the App's settings (More -> Your brush).

Note: Oral-B App functionality is limited if *Bluetooth*® is deactivated on your smartphone (for directions please refer to your smartphone's user manual).

- First pairing is indicated via the blue Smart Pressure Sensor light. Any App instructions will be displayed on your smartphone.
- First pairing needs to be done via Oral-B App. The toothbrush is automatically connected with your smartphone for future brush events.
- *Bluetooth* is enabled on your toothbrush by default. You may enable/ disable *Bluetooth* at any time by pressing the on/off button for 5 sec. Activation of *Bluetooth* is indicated via the light ring by a blue fade-in illumination. Deactivation is indicated by a blue fade-out illumination.
- Keep your smartphone in close proximity during brushing event. Make sure the smartphone is safely positioned in a dry place.


Note: Your smartphone must support *Bluetooth*® 4.2 (or higher) / *Bluetooth*® Smart in order to pair your toothbrush. Check [app.oralb.com](http://app.oralb.com) for compatible smartphones and country availability of the Oral-B App.

Caution: Follow your smartphone's user manual to ensure your phone/charger is designed for use in a bathroom environment.

## Cleaning recommendations

After brushing, rinse your brush head under running water with the handle switched on. Remove the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle as well as inside the brush head. Wipe all parts dry and let them air-dry before reassembling your toothbrush again.

Unplug the charger before cleaning. Clean the charger regularly with a damp cloth. Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.
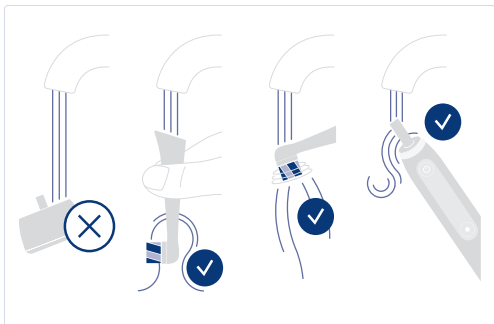
## Accessories (depending on model)

### Travel Case (d)

Store the clean and dry toothbrush/brush heads in the travel case (d).

The outer and inner case should be cleaned with a damp cloth only.

Note: Clean stains on travel case surface immediately. Store travel case in a clean and dry place.

### Brush head holder stand (e) / Brush head holder with cap (c)

You may store your cleaned brush heads on the pins of the brush head holders. Make sure brush heads are dried before

placing them in the holder and/or closing the cap. Brush head holders and cap are dishwasher-safe.

### Zipper Case / Magnetic Pouch

Store the clean and dry toothbrush/brush heads in the Zipper Case / Magnetic pouch.

The outer and inner parts can be cleaned under running water.

Store the clean and dry toothbrush along with the brush head holder (c) in the magnetic pouch. Make sure brush heads are dried before placing them into the holder and/or closing the cap. The outer and inner surfaces can be cleaned under running water. Store pouches in a clean and dry place.

Note: Clean stains on pouch surfaces immediately.

Subject to change without notice.


## Environmental Notice


Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

### Battery removal at the end of the product's useful life

Do not disassemble the product except when disposing of the battery. Run down battery before removing. When taking out the battery for disposal, use protective clothing (gloves, safety glasses, etc.) to omit skin contact whilst leaking batteries. Use caution not to short the positive (+) and negative (-) terminals. Insulate terminals (e.g. with tape) and put the battery in a plastic bag for disposal.

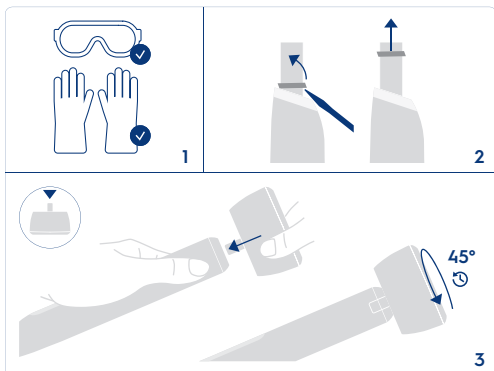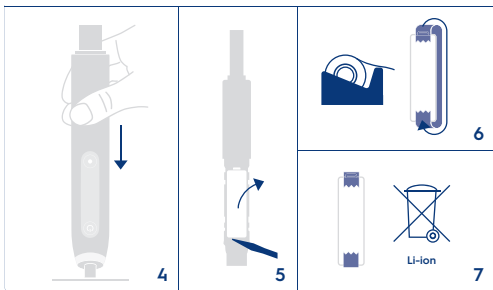Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

**Caution! Opening the handle will destroy the appliance and invalidate the warranty.**

## Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase.

Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different colour or an equivalent model.

This guarantee extends to every country where this appliance is supplied by Braun/ Oral-B or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun/Oral-B parts are not used.

This appliance is equipped with a *Bluetooth* Smart approved radio interface. Failure to establish a *Bluetooth* Smart connection to specific smartphones is not covered under the device guarantee, unless the malfunction of the *Bluetooth* interface is referable to Oral-B. The *Bluetooth* devices are warranted by their own manufacturers and not by Oral-B. Oral-B does not influence or make any recommendations to device manufacturers, and therefore Oral-B does not assume any responsibility for the number of compatible devices with our *Bluetooth* interface.

Oral-B reserves the right, without any prior notice, to make any technical modifications or changes to system implementation of device features, interface, and menu structure changes that are deemed necessary to ensure that Oral-B systems function reliably.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B/ Braun Customer Service Centre.

### **Warranty – Exclusions and Limitations**

This limited warranty applies only to the new products manufactured by or for The Procter & Gamble Company, its Affiliates or Subsidiaries («P&G») that can be identified by the Braun / Oral-B trademark, trade name, or logo affixed thereto. This limited warranty does not apply to any non-P&G product including hardware and software. P&G is not liable for any damage to or loss of any programs, data, or other information stored on any media contained within the product, or any non-P&G product or part not covered by this limited warranty. Recovery or reinstallation of programs, data or other information is not covered under this limited warranty. This limited warranty does not apply to (i) damage caused by accident, abuse, misuse, neglect, misapplication, or non-P&G product; (ii) damage caused by service performed by anyone other than Braun or authorized Braun service location; (iii) a product or part that has been modified without the written permission of P&G.

### **Information**

#### ***Bluetooth*® indicator**

Although all indicated functions on the *Bluetooth* device are supported, Oral-B does not ensure 100% reliability in connection, and feature operation consistency.

Operation performance and connection reliability are direct consequences of each individual *Bluetooth* device, software version, as well as the operating system of those *Bluetooth* devices and company security regulations implemented on the device.

Oral-B adheres to and strictly implements the *Bluetooth* standard by which *Bluetooth* devices can communicate and function in Oral-B toothbrushes.

However, if device manufacturers fail to implement this standard, *Bluetooth* compatibility and features will be further compromised and the user may experience functional and feature-related issues. Please note that the software on the *Bluetooth* device can significantly affect compatibility and operation.

### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable toothbrush damage is found to be attributed to the use of non-

Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable toothbrush at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

# Trouble Shooting

Problem	Possible reason	Remedy
<b>APP</b>		
Handle does not operate (properly) with Oral-B App.	<ol style="list-style-type: none"> <li>1. Oral-B App is closed.</li> <li>2. <i>Bluetooth</i>® is not activated on smartphone.</li> <li>3. <i>Bluetooth</i>® is disabled on handle.</li> <li>4. <i>Bluetooth</i>® connection to smartphone has been lost.</li> <li>5. Your smartphone does not support <i>Bluetooth</i>® 4.2 (or higher)/ <i>Bluetooth</i>® Smart.</li> <li>6. Outdated Oral-B App.</li> </ol>	<ol style="list-style-type: none"> <li>1. Launch Oral-B App.</li> <li>2. Activate <i>Bluetooth</i>® on smartphone (described in its user manual).</li> <li>3. Enable <i>Bluetooth</i>® by pressing on/off button for 5 seconds.</li> <li>4. Pair handle and smartphone again via App settings. Keep your smartphone in close proximity when using it with handle.</li> <li>5. Your smartphone must support <i>Bluetooth</i>® 4.2 (or higher)/ <i>Bluetooth</i>® Smart in order to pair handle. Check <a href="http://app.oralb.com">app.oralb.com</a> for compatible smartphones and country availability of the Oral-B App.</li> <li>6. Download latest Oral-B App.</li> </ol>
Reset to factory settings.	Original functionality wanted.	Use App to perform factory reset.
<b>TOOTHBRUSH</b>		
Toothbrush is not operating (during first use).	<ol style="list-style-type: none"> <li>1. Battery has a very low charge level.</li> <li>2. Operating outside standard temperature.</li> </ol>	<ol style="list-style-type: none"> <li>1. Charge at least for 30 min.</li> <li>2. Always operate in standard operating temperatures (&gt; 0°C and &lt; 40°C).</li> </ol>

<p>Toothbrush operates a short stuttering sound after 2 minutes only or does not operate any stuttering sound anymore.</p>	<p>Timer has been modified/ deactivated via App.</p>	<p>Use the Oral-B App to change timer settings or perform factory reset (see App Trouble Shootings).</p>
<p>Toothbrush does not charge.</p>	<ol style="list-style-type: none"> <li>1. Battery is at a very low charge level; it may take up to 30min until brushing mode display turns on.</li> <li>2. Ambient temperature for charging is out of valid range (<math>&lt; 0^{\circ}\text{C}</math> and <math>&gt; 40^{\circ}\text{C}</math>).</li> <li>3. Charger may not be (properly) connected.</li> </ol>	<ol style="list-style-type: none"> <li>1. Charge at least for 30min.</li> <li>2. Recommended ambient temperature for charging is <math>5^{\circ}\text{C}</math> to <math>35^{\circ}\text{C}</math>.</li> <li>3. Ensure charger is connected to mains voltage.</li> </ol>
<p>Rattling/loud noise during brushing or turning on.</p>	<ol style="list-style-type: none"> <li>1. Brush head may not be placed correctly onto the brush handle.</li> <li>2. Brush head may be worn.</li> </ol>	<ol style="list-style-type: none"> <li>1. Check proper fit of brush head on handle. Small gap remains in between.</li> <li>2. Replace brush head.</li> </ol>
<p>Mode button lights up in yellow colour.</p>	<p>Brush Replacement Reminder indicates that your brush might need to be exchanged.</p>	<p>Exchange Brush head and reset timer by holding mode button for 2 seconds.</p>

The rechargeable toothbrush contains a radio module which operates in the 2.4 to 2.48 GHz band with a maximum power of 1mW.

Hereby, Braun GmbH declares that the radio equipment types 3794, 3753, 3758, 3771, 3765, 3754, 3764, 3762, 3767 and 3776, are in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: [www.oralb.com/ce](http://www.oralb.com/ce).

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